

**Lab1**  
Universal Mixer

**One Mixer** - Endless applications...



# The adaptable in-container mixer designed for compounding

The Lab 1 universal mixer has been specifically designed for compounding. Fitted with a unique 3 dimensional inversion mixing arm which allows the user to obtain a perfect homogenized mix for both liquids & powder. The Lab 1 is also designed for “in container” mixing, any size\* any shape\*

Mixing the way that suits you, if traditional V-shell mixing is your preferred method then a V-shell container can be easily fitted in seconds ( no tools required ) giving you true flexibility to meet your customers needs.

Multiple vials, no problem the tube rotator attaches in seconds and allows gentle agitation of liquids by rotation.

## Features and Benefits

- Unique twin inversion motion for quick efficient mixing
- In container mixing
- Small quiet & compact
- Speed control from 0 to 60 rpm
- Interchangeable mixing arms for greater flexibility
- Adjustable angle from horizontal through 45 degrees
- Upgrade path for future attachments

## Technical specifications

Dimensions (base unit only): W x H x D ( in / mm )	(9 x 8 x 11 in) - (229 x 203 x 280 mm)
Weight (base unit only):	8kg
Speed range:	0 to 60 RPM
Electrical requirements:	110 - 240 V, 50/60 Hz
Maximum container size* Diameter x height:	(4 x 5 in ) - ( 101 x 127 mm )*
Warranty:	12 Months

# Test results for the Nymax Lab 1 universal mixer.

## Test conditions

Mixing device used: ( incl attachment type )	Lab 1 - Inversion arm
Substance mixed:	Pharmacy grade Lactose & Dyed corn starch
Container used:	750 ml jar ( material PET )
Mixing time:	90 seconds
Ambient condition:	Temp (celcius) 21 - Humidity 52%

Approximately 10g of dyed corn starch was placed into the mixing container with 60g of Lactose powder. ( fig -1 ) The container was placed into the mixer and mixed for 90 seconds at 60 rpm. ( fig - 2).



fig - 1

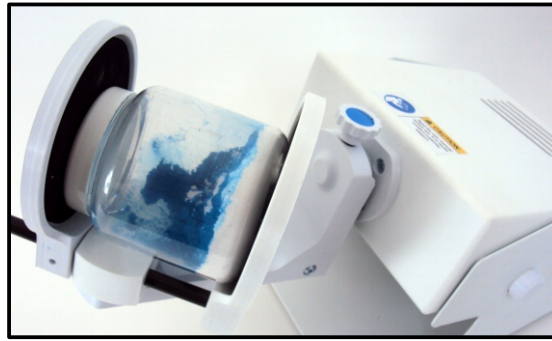


fig - 2

*Figure (1) ingredients in mixing container. Figure (2) Container being mixed for 90s at 60 rpm*

After mixing the powder was visually examined under a microscope at 50x magnification. ( fig - 4 )



fig - 3

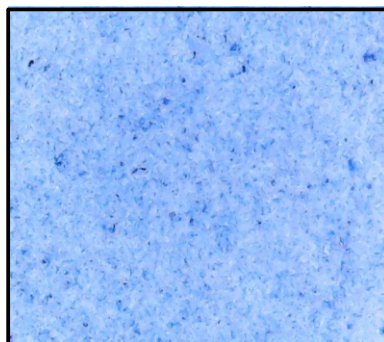


fig - 4

*Figure (3) Mixed ingredients shown in container. Figure (4) Mixed ingredients shown at 50 x magnification.*

## Conclusion

The Lab 1 mixer proved efficient at mixing powders resulting in a homogeneous mix.